## Z 09991

DE

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ab Seite 3

EN

### **Instruction manual**

starting on page 11

FR

### **Mode d'emploi**

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NL

### **Handleiding**

vanaf pagina 26

Z 09991\_V1\_12\_2015

## Inhalt

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## Content

Proper use _____	12
Safety notes _____	12
Product contents _____	15
Device overview _____	15
Before first use _____	15
Use _____	16
Technical data _____	17
Disposal _____	17

Dear customer,

We are delighted that you have decided to purchase this **tooth polisher**.

Before using the device for the first time, please carefully read through the instructions and store them for future reference. These instructions are to accompany the device when passed on to others.

The manufacturer and importer assume no liability in the event the data in these instructions have not been observed!

If you have questions about the device or spare parts/accessories, contact customer service via our website:

**[www.service-shopping.de](http://www.service-shopping.de)**

## Explanation of symbols

Safety instructions: Please carefully read through and obey the safety notes in order to avoid injury to persons and damage to property.

## Proper use

- This device is for polishing teeth. Do not use it on any other part of your body.
- The device is designed for private, not for commercial use.
- Only use the device as described in this instruction manual. Any other use is deemed improper. Store the instructions in a safe place.
- Defects arising due to improper handling, damage or attempts at repair are excluded from the warranty. This also applies to normal wear and tear.

## Safety notes

### Risks of injury

- Keep children and animals away from the packaging material. Danger of suffocation exists!
- This device is not suitable for use by persons (including children) with restricted sensory or intellectual abilities or with a lack of experience and/or knowledge unless they are supervised by a person responsible for their safety or have received instructions from that person as to how to use the device.
- Children must be supervised to ensure that they do not play with the device.
- Do not use the device if it is damaged. Each time before you use the device, check to make sure that the device and the adapters are in good condition. In the event of damage or wear, injuries may result.
- Polish all your teeth evenly and do not forcefully push the adapter against your teeth.
- Immediately stop using the device if it is unpleasant or painful when applied.

- To prevent accidents, keep the rotating adapter away from your hair, textiles etc.
- Make sure that the rotating adapter of the device only comes in contact with the teeth. Do not touch it with your fingers!

### **Fire hazard**

- Do not use the device if easily flammable gases are present in the air.

### **Proper handling of batteries**

- Batteries may be life threatening if swallowed. Therefore keep the batteries and device out of the reach of children and animals. If a battery has been swallowed, you must seek immediate medical attention.
- Before inserting the batteries, check whether the contacts in the battery compartment and on the batteries are clean and, if necessary, clean them.
- Only use the battery type specified in the technical data (see section “Technical data”).
- When replacing the batteries, make sure that the polarity (+/–) is correct.
- Remove the batteries from the battery compartment if they are empty or if you are not using the device for an extended period. This way, you can avoid damages attributed to leaking.
- If the batteries have leaked, avoid contact of the skin, eyes and mucous membranes with the battery acid. In the event of contact with the battery acid, rinse the affected areas instantly with plenty of clean water and immediately consult a physician.

- Immediately remove any leaking batteries from the battery compartment. Clean the contacts before you insert new batteries.
- The batteries may not be charged or reactivated with other means and they may not be disassembled, thrown in a fire, submerged in a liquid or short-circuited.


## Handling the device properly

- Check the device for damages prior to each use. If the device or the adapters exhibit visible damages, the device may not be used.
- Never attempt to repair an electrical device yourself, rather contact customer service centre or locate a qualified workshop to guarantee your own safety.
- Only use the device if the adapter has been correctly mounted!
- Only use the adapters from the manufacturer. Only use the adapters with this device.
- Do not allow the device to fall and do not subject it to strong percussions.
- Never submerge the device in water or other liquids.
- Do not use any corrosive or abrasive cleaning agents to clean the device. They could damage the surface. Clean the adapters with warm water and, when needed, wipe the device off with a clean damp cloth. Then dry it all off.
- Keep the device and the accessories away from sources of heat.
- Do not subject the device to extreme temperatures, severe temperature fluctuations, direct sunlight or moisture.
- Store the device and accessories in a dry place out of the reach of children and animals.

## Product contents

- 1 x main unit
- 1 x polishing adapter
- 1 x storage bag
- 2x 1.5 V type AA batteries

## Device overview

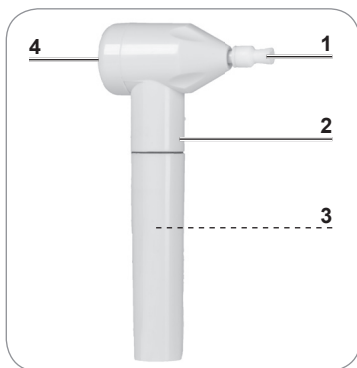- 1 Adapter
- 2 Housing
- 3 Battery compartment (in the housing)
- 4 On/Off switch

## Before first use

- Unpack all parts and make sure that the product contents are complete (see chapter “Product contents”) and check for any damages from transport. If applicable, remove any stickers and foils. If you identify damages to the components, do not use them (!), instead contact our customer service centre.
- Insert batteries (see chapter “Inserting batteries”).

## Inserting batteries

### ATTENTION!

- Always make sure that the device has been switched off before inserting/changing batteries.
- 
1. Turn the housing counterclockwise.
  2. Pull the housing off.
  3. Insert two 1.5 V batteries (type AA) in the battery compartment. Make sure the polarity (+/-) is correct.
  4. Put the housing back on the battery compartment. Make sure that the catches are located in the designated recesses.
  5. Turn the housing clockwise to tighten it.
- Change the batteries as soon as you notice a drop in power.

## Use

### ATTENTION!

- Only use the device on your teeth. Do not use it on your skin!
- Polish all your teeth evenly and do not forcefully push the adapter against your teeth.
- Immediately stop using the device if it is unpleasant or painful when applied.

## Attaching/removing adapters

- Attaching: Mount the adapter on the pin (1) with the wider end forward.
- Removal: Simply pull the adapter off of the pin.

## Using the device

We recommend polishing your teeth no more than twice a month in order to protect them.

1. Make sure that your teeth are clean.

2. Turn the device on by pushing the On/Off switch (4).
3. Place the rotating adapter up against a tooth and slowly move it on the tooth without applying pressure until it shines.
4. Once you have polished all teeth, Turn the device off by pushing the On/Off switch again.

## Technical data

Model number: MF11/MC0085

Voltage supply: 3V DC (2x 1.5V batteries, type AA)

## Disposal

The packaging material is recyclable. Dispose of the packaging in an environmentally-friendly manner and make it available to the collection service for recyclable materials.

Dispose of the device in an environmentally-friendly manner. It does not belong in the household garbage. Dispose of it at a recycling centre for old electric and electronic equipment. For more information, please contact the administration in your community.

Batteries and rechargeable batteries must be removed and disposed of separately before disposing of the device. To protect the environment, batteries and rechargeable batteries may not be disposed of together with normal household waste. Instead they must be turned in at appropriate collection points. Please also observe applicable statutory provisions regulating the disposal of batteries.

### Customer service/importer:

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Am Heisterbusch 1, 19258 Gallin, Germany

Tel.: +49 38851 314650 \*)Calls subject to a charge. Costs vary depending on the provider.

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All rights reserved.
DS Produkte GmbH