




Model Felicitas / BM 702 Instruction Manual



Directory

- 1. Conditions of use.....page 2
- 2. Commissioning.....page 2
- 3. Weighing function.....page 3
- 4. Technical Data.....page 3
- 5. Warranty.....page 3

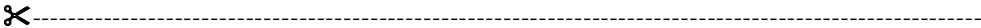

Guarantee card – 3 years warranty for the Mechanical Person Scales

<ul style="list-style-type: none"> • Sender <hr/> <hr/> <hr/>	<ul style="list-style-type: none"> • Reason for complaint <hr/> <hr/> <hr/> <hr/>
--	--

• Date of purchase

Please return this guarantee card together with the receipt to your retailer.

Thank you for purchasing this high-quality product made by ADE. You have acquired a set of scales that meet the highest technical standards and are of proven quality.

Please read these operating instructions carefully before using the product for the first time and keep them in a safe place where you can find them whenever you need the information contained in the instructions.

2. Conditions of use

If possible, the weighing should be done undressed. Please keep the weighing conditions constant always at the same daytime.

By nature, a body is subject of fluctuations (for example loss of water after taking exercise or after a sauna visit, after a meal, etc.). To enable comparable readings it is important to take the weighing, if possible, always under the same conditions.

ADE takes no responsibility for any damages or losses caused by the use of this scales or for Third Party claims. This product is exclusively intended for private use. This product is not designed for commercial or professional use in hospitals or other medical facilities.

2. Commissioning

Before first time use of the scales, please check at the beginning carefully if there are any damages at the packing or at the device. Not on any account use the device if it is damage. In this case, please immediately contact us. Our staff members willingly will help you.

Safety guidelines:

The appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.



Place the scales on a plane and firm surface (tiles, hardwood flooring, etc.) to avoid a wrong measurement on carpeted floor.



Always stay in the centre of the scales to prevent the scales tipping to the side.



Attention - Hazard of slips! Do not stand on the scales with wet feet or on a slippery surface of the weighing platform.



Please care that no heavy object may fall on the scales. The scales might become damage.

Clean the scales only by using a slightly damp cloth. Do not use solvents or abrasive cleaning agents. Do not submerge the scales under water.

3. Weighing function

1. On the lower edge of the scales there is a small rotary knob for precisely setting the zero point. Please make sure that the hand of the scale reading at no load shows exactly on the zero point.
2. Step on the centre of the scales. Remain stationary on the weighing platform and do not hold you tight. After a few seconds, your weight will be indicated.

4. Technical Data

Capacity x graduation: max. 160 kg x 1000 g
Dimensions: 275 mm x 345 mm x 113 mm

Subject to technical modifications

5. Warranty

ADE warrants for a period of 3 years from the date of purchase that it will remedy any defects due to faulty material or workmanship free of charge by repair or replacement. When purchasing, please have the warranty coupon signed and stamped by the dealer. When making a guarantee claim, please return the weighing scale together with the guarantee card including the reason of complaint to your dealer.

