PUREFX

HIGH EFFICIENT MASSAGE DEVICE

USER MANUAL

ACTIVE  HERO

Summary

Get started with PURE FX	3
Packaging content.....	5
Description of the PURE FX	6
Getting started.....	7
Instructions for use.....	9
Different massage heads.....	10
Different vibration speeds.....	10
Areas of use of the PURE FX	11
Advice on care, safety and tips for use.....	12
Technical features.....	14
Disposal.....	14

Massage guns are now an integral form of sports reinforcement exercises and are used by sports professionals, weight lifters and endurance athletes.

Get started with PURE FX

Physiotherapists, chiropractors and sports coaches are turning to vibration therapy to help athletes improve their performance and recover from their injuries. **PURE FX** is your new fitness ally.

Use **PURE FX** before your workouts to increase your performance or after your workouts to speed up recovery. The vibration technology of the **PURE FX** penetrates deeply into your muscles and fascias.

The high frequency pulsations of the **PURE FX** can relax contracted muscles, increase the blood supply to the affected area and reduce inflammation.

The **PURE FX** is easy to carry and its lithium-ion battery provides use for almost 3 hours. Different massage speeds and pressure intensities allow you to select a power adapted to each of your needs.

Give your training a new dimension and discover the power of PURE FX!

Go through the exercises described in this booklet to find out how to get started with your **PURE FX**.

PUREFX

HIGH EFFICIENT MASSAGE DEVICE

Packaging content:

- PURE FX – High efficient Massage Device
- charger
- transport case
- 4 x massage heads
- instructions for use

Please read this user manual carefully, keep it for future use, make it available to other users and follow its instructions. Follow the steps in the **Getting started** section on page 7 to mount and load your PURE FX.

Description of the PURE FX

LED display of the speed level

The LED displays the current speed levels.

Power button

Press the power button to turn on the **PURE FX** and set it to speed level 1.

Press the power button again to increase the speed level.

Pressing the activation key five times turns the **PURE FX** off again.

LED display for the current charge level

Getting started

Fully charge the battery

- Fully charge the **PURE FX** before first use (up to six hours)
- To do this, connect the power supply to the **PURE FX** charging socket
- Plug the adapter into a power outlet
- When the battery is fully charged, unplug the power supply

Attaching a massage head

- Once your **PURE FX** is fully charged, you can attach a massage nozzle
- Choose the desired massage nozzle
- Push the nozzle firmly inside to secure it
- The spinal nozzle has a small nodule. Place the base of the nozzle so that it is facing the notch on the nozzle of the **PURE FX**

Also see the different **massage heads and vibration speeds** on page 10.

PURE FX

Instructions for use

Read the instructions before using the **PURE FX**. Keep this manual for future reference.

Charging

- Completely charge the **PURE FX** before its first use (up to 6 hours)
- After the first charge, the **PURE FX** will require no more than around 3 hours charging time
- The battery is fully charged when the LED on the power supply goes from red to green
- Maximum usage time extends up to around 4 hours (depending on the speed used)

Using **PURE FX**

- Push the ON / OFF switch to ON to turn on the **PURE FX**
- Press the activation switch to start up the **PURE FX** at speed level 1. Select one of 4 levels of vibration speed: Press the activation switch to increase the vibration speed
- Move the head of the massage nozzle along a muscle group at the desired speed level. Exercise different levels of pressure according to your needs.

Removing and replacing a massage head

- Push the ON / OFF switch to OFF to turn off the **PURE FX**
- Grasp the attached nozzle and pull firmly to remove it
- Insert the base of another nozzle into the device
- For the spinal nozzle, be sure to place the small nodule of the nozzle facing the notch on the nozzle of the device
- Push the nozzle firmly inside to secure it

TIPS FOR USING THE **PURE FX**

- Move the **PURE FX** over the desired muscle groups for between 15 seconds and two minutes
- Adapt the pressure applied according to your needs, but never use excessive force
- Never use the massage head for more than 3 minutes on the same area!
- Prolonged use could cause tissue damage

Different massage heads

Four interchangeable massage nozzles are included with your **PURE FX**. These massage nozzles are designed to treat different areas of your body and meet specific needs for therapeutic massage.

Ball

For recovery and reduction of body aches after training. Suitable for all muscle groups.

Flat

For muscles and areas that need special attention. Simulates a deep massage.

Cone

Specifically designed for problem areas and the treatment of « trigger points ».



Spinal

Ideal for back muscles situated along the spinal column.


Areas of use for the PURE FX

For all Muscle groups

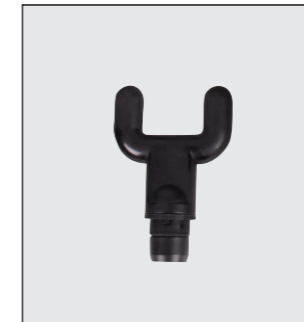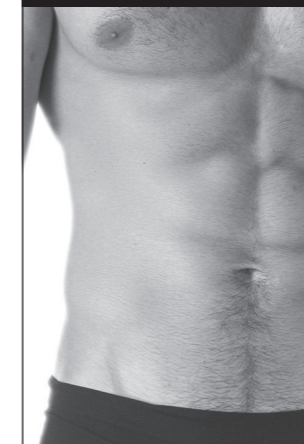For back muscles

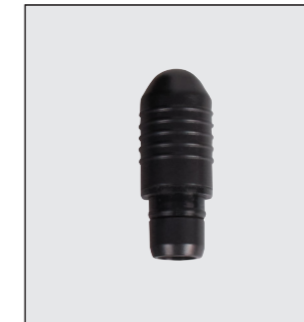For punctual treatment

For all parts of the body

Different vibration speeds

The 4 different massage speeds make it possible to select the right intensity for a specific muscle group or trigger point. Each speed offers a different vibration frequency (Hz) and a different number of revolutions (Rpm).

The number of revolutions can vary between 1,800 and 3,200 Rpm.

Advice on care, safety and tips for use

Maintenance, cleaning and storage

- Use a soft cloth to wipe off dirt or moisture on the device
- When you want to store or transport it, remove the battery and store it in the included case

Safety notes

- Please read the following warnings carefully to minimize the risk of burns, fire, electric shock or injury
- **Seek medical advice before using this product. Do not use it if you are injured**
Exclusively for adults. Keep out of the reach of children and pets
- **The use of this device is not a substitute for a medical consultation or treatment.**
In case of any pain or illness, consult your doctor immediately!
- Use the **PURE FX** exclusively according to the instructions provided in this manual
- Keep away from liquids and heat sources
- Never let the **PURE FX** run unattended
- Unplug the appliance from the mains after charging or before use
- Do not remove any screws or attempt to disassemble the device
- Do not exceed one hour of consecutive operation. After one hour of operation, leave the unit off for 30 minutes before using it again

For the sake of your health, it is not recommended that you use the massage gun in the following cases:

- Presence of electrical implants (such as a pacemaker, for example)
- Presence of metal implants
- Whilst wearing an insulin pump
- In case of fever
- In case of known or acute heart rhythm disorders and other disorders of cardiac excitability and conduction
- In case of seizures (eg epilepsy)
- Pregnant women
- People suffering from cancer
- Following operations where increased muscle contractions may interfere with recovery
- Never use the device near the heart. Stimulation electrodes should never be placed in front of the ribcage (the ribs and sternum of which form the visible part), and especially not on the two large pectoral muscles. Such use may increase the risk of ventricular fibrillation and lead to cardiac arrest
- On the skull, in the mouth, throat or larynx
- On the neck / carotid
- On the genitals
- On skin with severe or chronic lesions or inflammations, eg. in case of inflammation, redness, rash, burns, bruises, swelling, painful or not, and in case of open wound or healing, on postoperative scars during cicatrisation
- In environments with high humidity, eg. in a bathroom, in the bath or in the shower
- Do not use after consuming alcohol

Important precautions for use

- Only use **PURE FX** on soft tissues
- Do not use **PURE FX** on your head or other bone areas
- Hematomas may appear, regardless of the chosen rate of vibrations
- Control the areas you treat frequently and adapt your use at the first signs of pain or discomfort
- Use only on clothes or clean dry skin
- We recommend that you use the **PURE FX** over your clothes (eg. tight fitting sportswear, leggings, etc.)
- Keep your hair, fingers, and other body parts away from the base of the massage head to prevent bruising

WARNING!

Use this device exclusively for massage:

- On people
- For the purpose for which it was developed and in the manner indicated in this instruction manual
- Incorrect use can present risks
- External use
- With the original accessories supplied and re-orderable, otherwise the right to the guarantee will not apply

Technical data

Article number:	21402 / HERO02
Model:	D06
Dimensions:	Approximately 17 cm x 23 cm x 7 cm
Battery:	Lithium-Ion
Battery lifetime:	Up to 4 hours (depending on the speed used)
Battery charging time:	Around 3 hours
Power supply:	220 - 240 V ~ 55 Hz
Power:	1500 mAh 24 V
Speed:	4 speed levels up to 3200 Rpm
Noise level:	Approximately 56 DB at maximum Speed
Weight:	Approximately 1.2 kg (battery incl.)
protection class:	II
Made in PRC	

Disposal

Packaging materials are recyclable. Dispose of the packaging in an environmentally friendly way and take it to a recycling centre.

Dispose of the device in an environmentally friendly manner. It should not be thrown away with household rubbish. Dispose of this unit in a recycling centre for used electrical and electronic equipment. For more information, consult the local guidelines of your borough.

Broszio & Co.
GmbH & Co. KG
Ottensener Str. 1-5
D-22525 Hamburg
Tel.: 040 - 5472 92 0
E-Mail: info@broszio.com
www.broszio.com

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